

# Installation Instructions

## Stacking Racks

Part #AVA50500, #AVA50505, #AVA50510

### Description:

The Apollo Dental Products Stacking Rack provides a simple method of saving crucial floor space by placing the vacuum system directly above the air compressor. The Stacking Rack is a high strength, yet lightweight structure that is easily assembled with the use of a soft metal or rubber hammer.

### Assembly Instructions:

1. Remove the top plate assembly, leg assemblies, spring clips and 3/8" tubing from the shipping enclosure.
2. Place the top plate assembly on the floor so the corner ears are facing up.
3. Position two spring clips in the indentations of a corner ear and slide a leg assembly over it (Figure 1).
4. With a soft metal or rubber hammer, tap the leg assembly until it is fully engaged onto the corner ear.
5. Repeat steps 3 and 4 for the remaining three leg assemblies.
6. Attach the 3/8" safety drain tubing to the adapter located on the underside of the top plate assembly.
7. Position the rack over the compressor.
8. Adjust the leveling feet so that the rack rests firmly on all four feet and has a slight slope towards the rear.

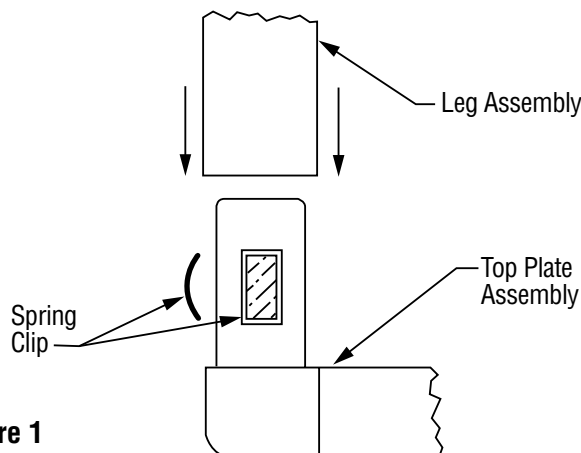
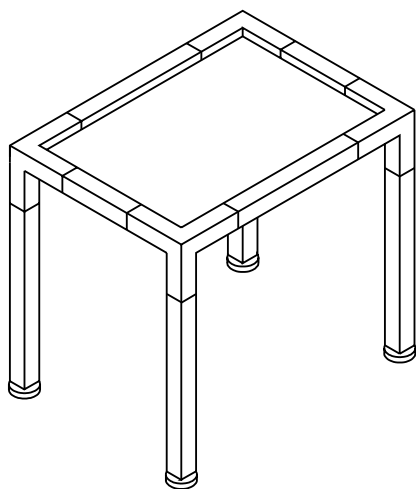


Figure 1

### Dimensional Specifications:

Model Number	Width (In.)	Depth (In.)	Height (In.)	Weight
AVA50500	26	28	37	26.5 lbs.
AVA50505	40	28	37	32.5 lbs.
AVA50510	44	30	40	36 lbs.

### Warranty Information: 90 Days

All ADP accessories are thoroughly inspected and tested in accordance with rigid specifications and standards. Our accessories are guaranteed against any defective material and workmanship from the date of shipment; provided that the installation, operation, and maintenance is done in accordance with ADP procedures as outlined in our Installation and Maintenance Guides. Warranty cards must be returned to ADP within ten days of installation to effect warranty. No other warranties or guarantees, expressed or implied are made.

ADP's obligation under the warranty is to provide parts for the repair or, at its option, to provide the replacement product (excluding labor). All special, incidental and/or consequential damages are excluded. We will not issue credit for accessories without first attempting to correct the problem in the field. Written notice of breach of warranty must be given to ADP within the warranty period. The warranty does not cover damage resulting from improper installation or maintenance, accident or misuse. The warranty does not cover damage resulting from the use of cleaning, disinfecting or sterilizing chemicals and processes. The warranty does not cover vacuum failures due to hard water deposits. Failure to follow instructions provided in ADP's Installation and Maintenance Guides may void the warranty.